


Expand your clinical solutions with the

TUFF UNICON

Conical Connection Implant System

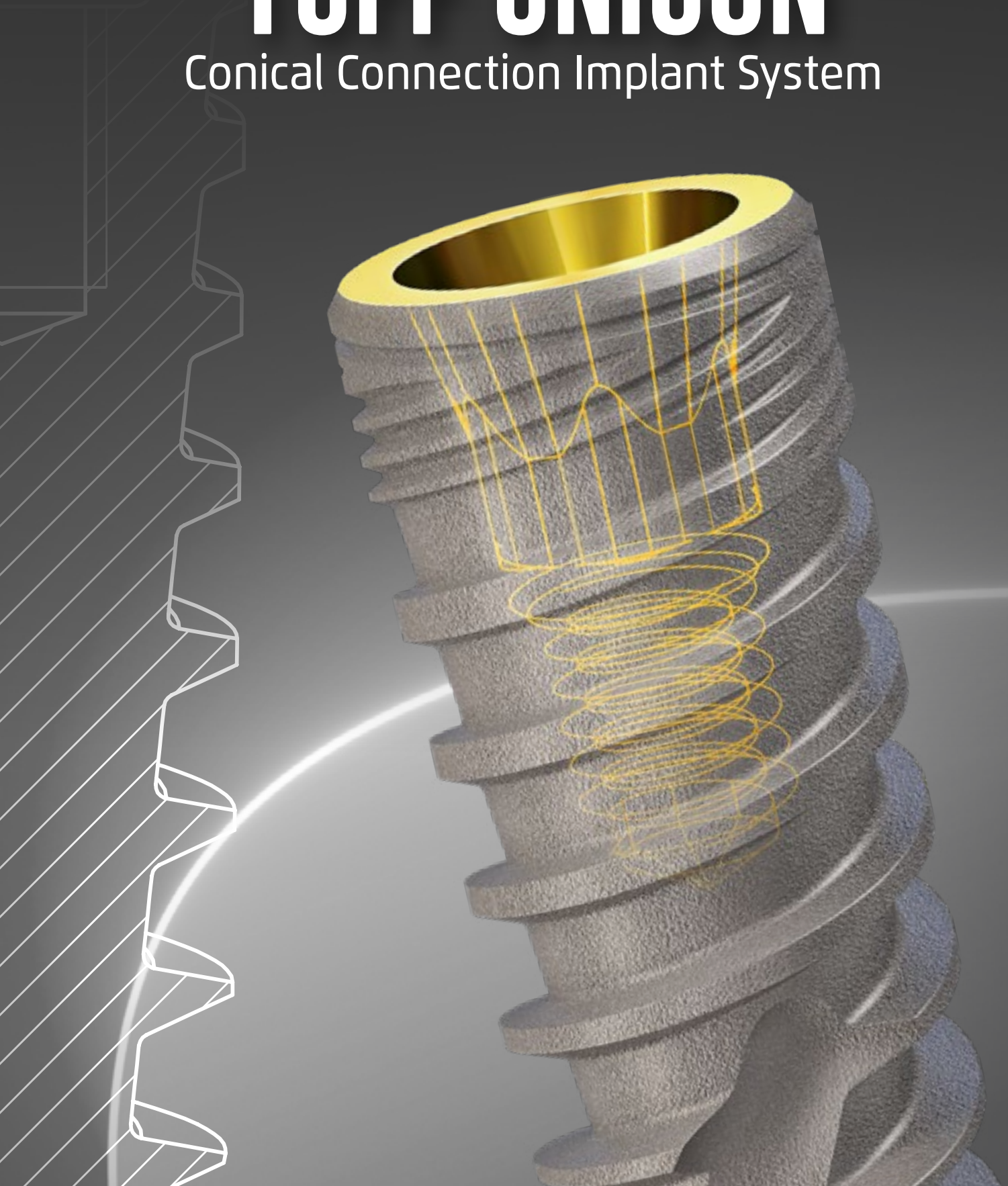

Expand your clinical solutions and provide your patients with the best in dental implant technology

FOR THE BEST STABLE AND AESTHETIC RESULTS

Our TUFF Unicon implant has all the benefits and features of our excellent TUFF implant with a new conical connection



CONICAL CONNECTION

- High fatigue strength due to load distribution
- Tight implant-abutment fit & Minimal micro-movements
- Large platform switching

CORONAL MICRO THREADS

- Decrease the load transfer to the crestal cortical bone
- Increased surface area


LARGE STEP DOUBLE THREAD DESIGN WITH TAPERED CORE BODY

- Efficient and fast insertion
- Maximum dissipation of stress

THE LOWER V-SHAPE THREAD ZONE

Fast & smooth insertion with excellent retention for immediate loading

Conical Prosthetic Solutions
**ENHANCING SOFT TISSUE MANAGEMENT,
AESTHETIC EXCELLENCE AND LONG-TERM STABILITY**

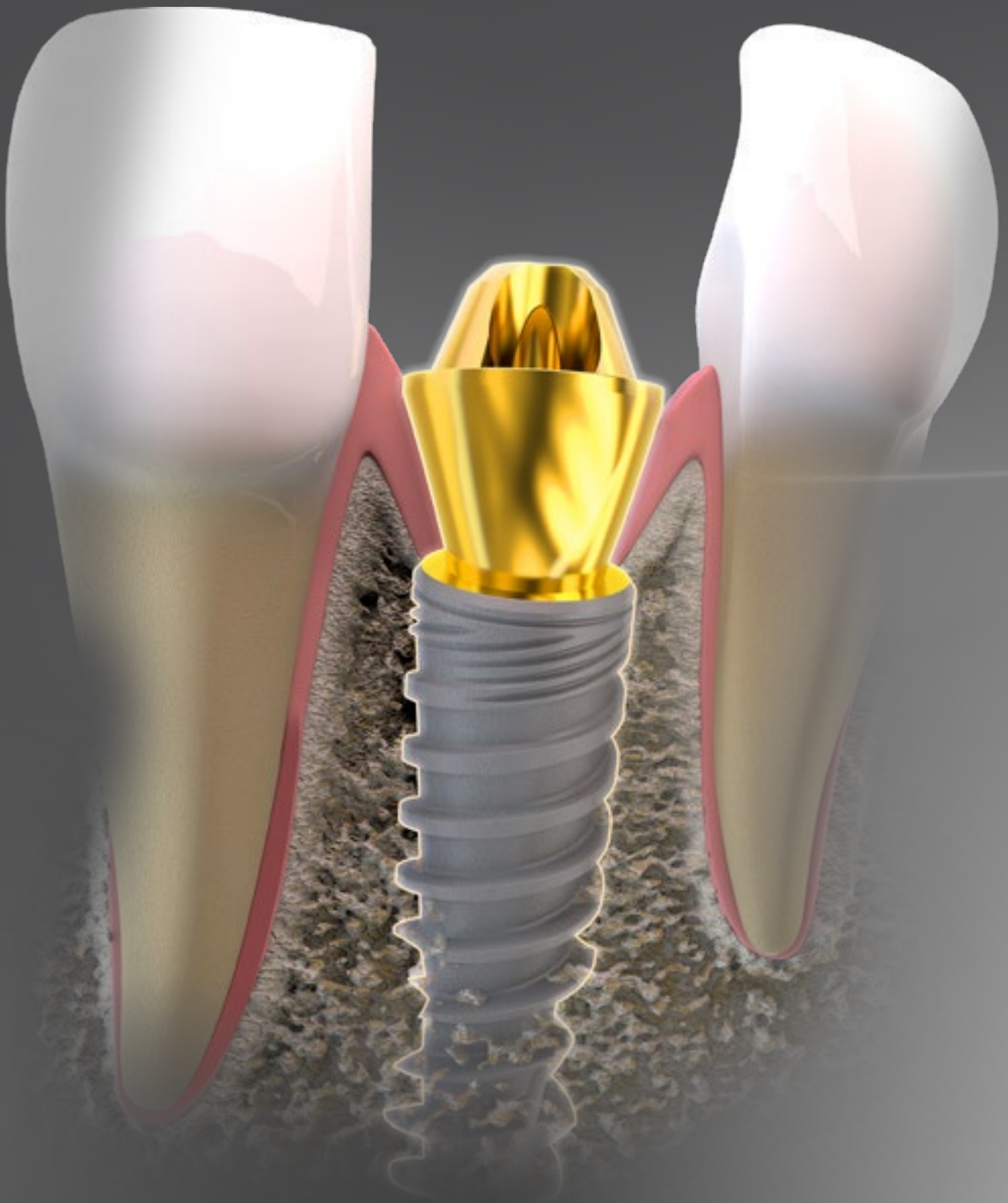

CONICAL CONNECTION PLATFORM




DESIGN FEATURE

- Small angle 24° Internal Cone



CLINICAL BENEFITS

- High fatigue strength due to load distribution
- Stability of prosthetics for minimal micro movements
- Optimal sealing to minimize bacteria/food residue leak
- Emergence profile angle that allows soft tissue space

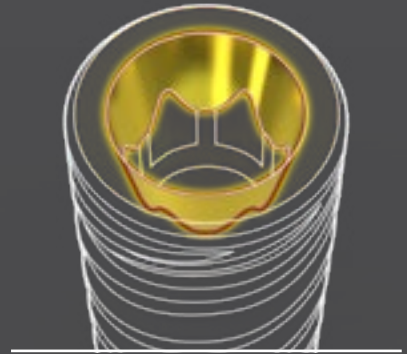


DESIGN FEATURE

- 2.26 mm Internal Hexagon



CLINICAL BENEFITS

- 6 positions simple/functional for prosthetics orientation

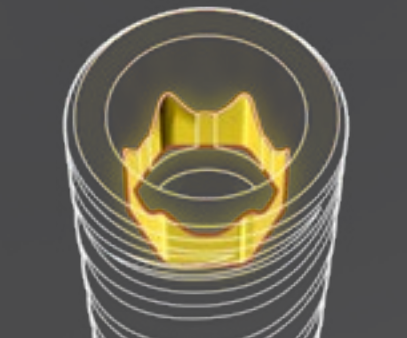

DESIGN FEATURE

- Flat coronal upper surface




CLINICAL BENEFITS

- Large Platform Shifting/Switching
- Flat top surface for accurate sealing of healing abutments
- Allowing space for soft tissue to grow



CONICAL CONNECTION PLATFORM


Strong and Stable Implant-abutment connection with a small angle for strength and optimal tissue engagement

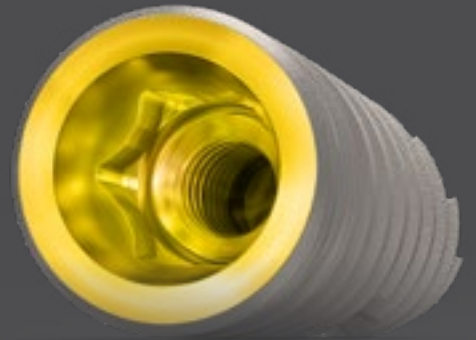

ONE PLATFORM FOR ALL IMPLANT DIAMETERS

- Same Hexagon
- Same angle
- Same internal thread

*For $\varnothing 3.25$ Implant diameter smaller interface diameter for the same prosthetics

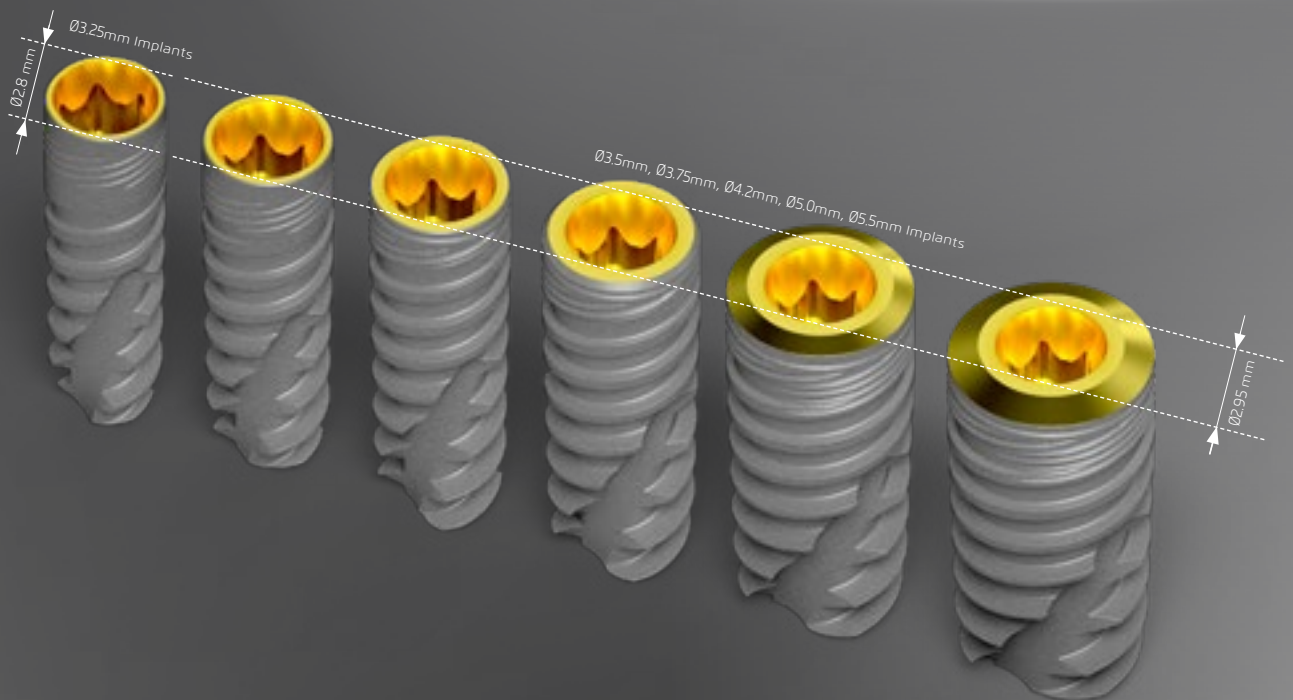

Conical Prosthetic Solutions

CONICAL CONNECTION SOLUTION FOR SUPERIOR STABILITY AND AESTHETICS

Noris Medical offers a comprehensive solution utilizing the conical connection platform.

This intuitive, state-of-the-art conical connection solution is suitable for most clinical indications, ranging from single tooth replacements to fully edentulous cases.

The Complete system perfectly integrates biological principles with an aesthetic approach.

EVERYTHING YOU NEED IN ONE PLATFORM

from Ø3.25mm implant to Ø5.5mm
one platform, one prosthetic solution, one kit.

As an all-inclusive conical connection solution, the TUFF Unicon concept seamlessly integrates implants, prosthetics, the Noris Medical surgical tool kit and robust digital workflows.

This harmony between the predictability and performance of our conical connection implant system and our advanced prosthetic solutions ensure superior results.

CONICAL CONNECTION IMPLANT SYSTEM



PROSTHETIC PRODUCT LINE



IMPLANT DRIVERS



PREMIUM SURGICAL SET




Conical Prosthetic Solutions

ENHANCING SOFT TISSUE MANAGEMENT, AESTHETIC EXCELLENCE AND LONG-TERM STABILITY



MULTI-UNIT ANATOMIC/CONCAVE SCREW RETAINED ABUTMENTS

Professionally designed abutments for optimal tissue management, aesthetics, and patient comfort for amazing aesthetic results.

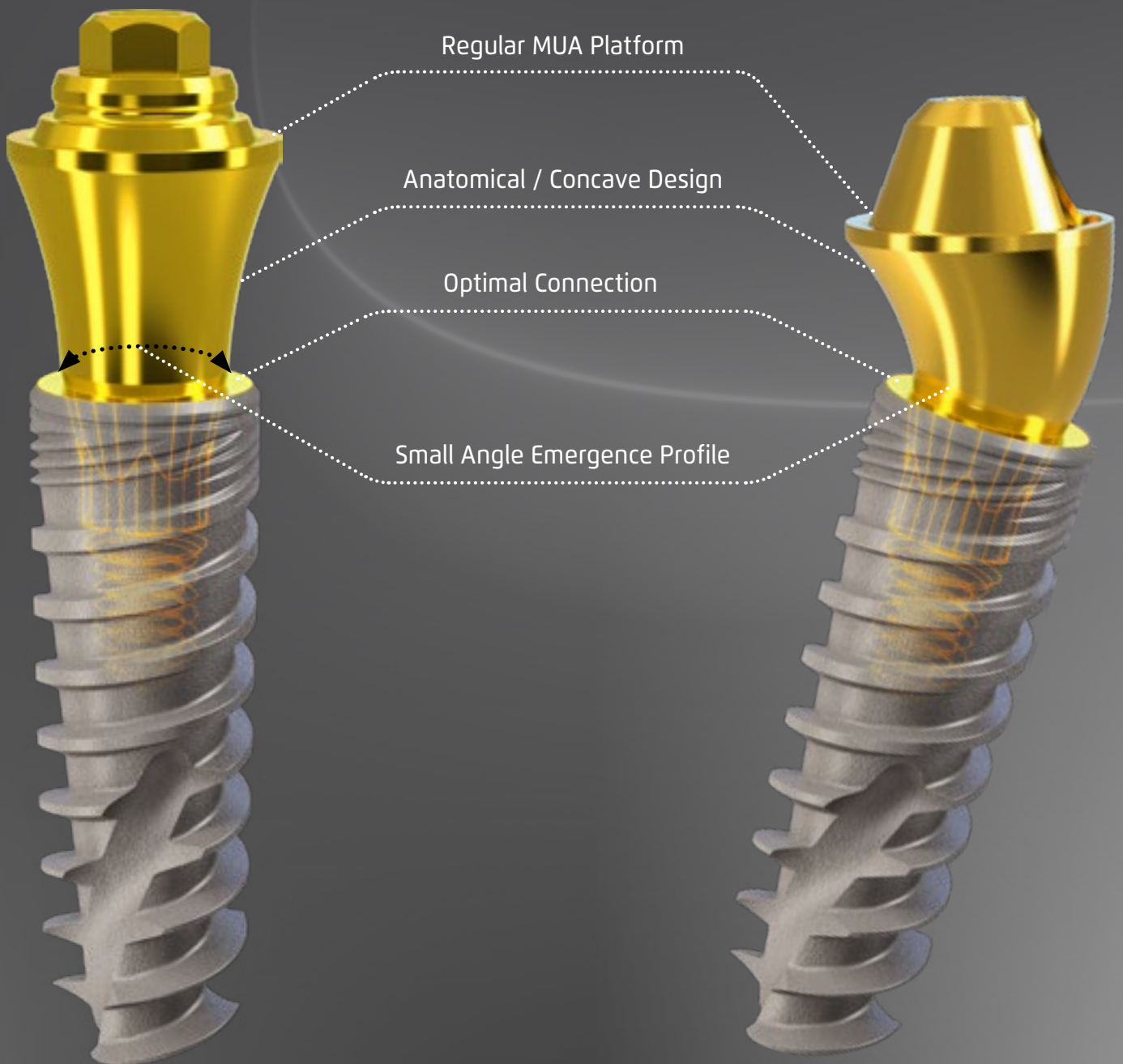

TUFF UNICON PROSTHETIC PRODUCT LINE

Enhancing Soft Tissue Management, Aesthetic Excellence and Long-Term Stability

HEALING CAPS

- Wide range of gingival heights and diameters to fit all prosthetic products
- Diameters: Ø4.3, Ø4.8, Ø5.2
- Allow healing of the surrounding soft tissue
- Support the healing process
- Prevent gum overgrowth
- Reduce risk of infection



IMPRESSION TRANSFERS

- Closed Tray and Open Tray
- Available in short, high, narrow and wide sizes



ANALOGS

- For cast model
- Gold Anodize to represent the Conical Connection platform
- Representing: Ø3.25, Ø3.5, Ø3.75, Ø4.2, Ø5.0, Ø5.5 implants



TUFF UNICON PROSTHETIC PRODUCT LINE

Enhancing Soft Tissue Management, Aesthetic Excellence and Long-Term Stability

TEMPORARY ABUTMENTS

- Gold Anodize to represent the Conical Connection platform
- For bridge or crown restoration
- Deep retention grooves



SCREW RETAINED ABUTMENTS

- Gold Anodize to represent the Conical Connection platform
- The same MUA platform as Internal Hex MUA
- Anatomical emergence profile
- Optimal esthetic results



CEMENT RETAINED ABUTMENTS


- Gold Anodize to represent the Conical Connection platform



TOOLS

IMPLANT DRIVERS

- Self-holding
- Available in various lengths

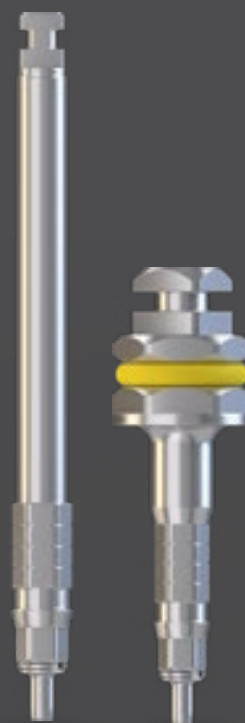

PREMIUM SURGICAL SET



We can make you
SMILE